


Frequency Comb

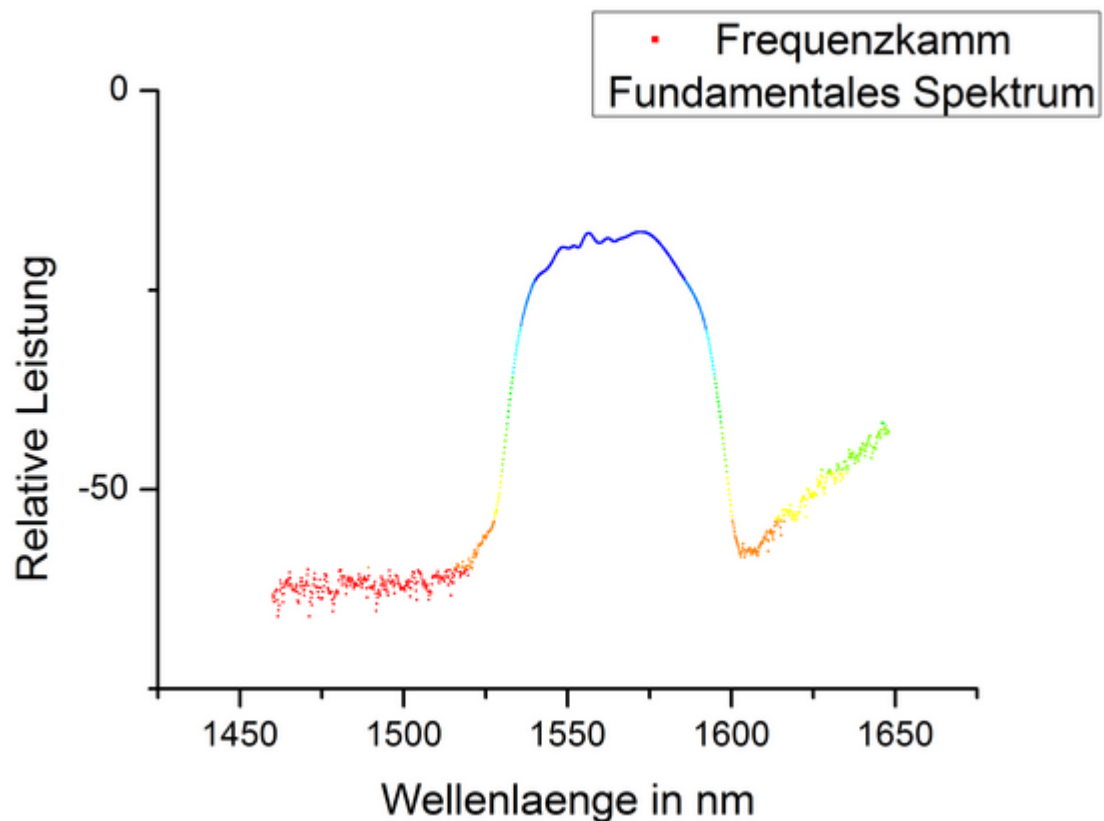


Setup


Typical Values

- Laser diode current: 9.0
- AC: 950, DC: 850 (in Fiber laser software)
- Powers:
 - after fiber laser: ~ 200 mW @ 1560 nm
 - after SHG:
 - ~ 3 W @ 780 nm
 - ~ 95 mW @ 1550 nm
 - direct after filter:
 - ~ 32 mW @ 780 nm
 - direct infront PCF:
 - ~ 300 mW @ 780 nm without filter
 - ~ 20 mW @ 780 nm with filter
 - ~ 8 mW @ 1550 nm without filter
 - after PCF:
 - ~ 8.5 mW @ 780 nm without filter (red light should be visible)

- $\sim 1.8 \text{ mW}$ @ 780 nm with filter
- $\sim 250 \mu\text{W}$ @ 1550 nm without filter
- Clock Laser / DUT $\sim 5 \text{ mW}$; FWHM $\ll 100 \text{ MHz}$
- (right picture is an old measurement and we don't know where in the comb)
- CEO frequency has a linewidth of $\sim 500 \text{ kHz}$
 - yellow = good lock
 - white = bad lock, because the platform and width

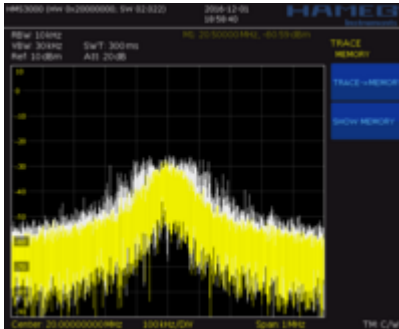

Electronic part

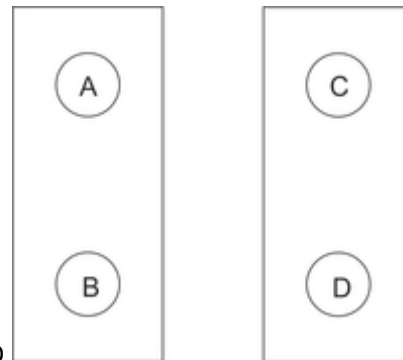

- Die Reprate kommt vom Synthesizer: $10 \cdot f_{\text{rep}}$
- $10 \cdot f_{\text{rep}}$ geht in die Reprate-Box [A]
- Photodioden Signal des Frequenzkammes von "RF OUT"
- Dort wird [A] und [B] gemixt
- Rapraten PID steuert Piezo im Resonator (Länge: 1.5 m $\sim$ FSR = 100 MHz)
- Reprate ist fouriertransformiert der zeitliche Abstand zwischen zwei Pulsen

Regular problems

- Comb has to be optimized almost every day: Most of the times just tuning potis @ SHG & PCF is enough (see picture)



- Low power in general: Maybe pump diodes should be exchanged (15.09.2016: 145 mW Output @ 9.5)
- Fiberlaser software does not start: Maybe wrong COM Port: Change @ Windows Device Manager
- FiberLaser software starts, but Fiber Comb does not: Maybe NI card not attached properly
- You can't lock CEO and Reprate, than restart the hole lasersysteme [probable all ~1-2 months].


Turn off and on



- Turn off the frequency comb:
 1. Push Amplifier button (LED off)
 2. Push Laser button (LED off)
 3. turning the upper key (position 0)
 4. turning the lower key (position 0)
- Turn on:
 1. exactly the other way around

1064 Port


Data under ResLab: D:\Dokumente\Daten\Frequenzkamm\Spektren 161115 1064 nm Port.

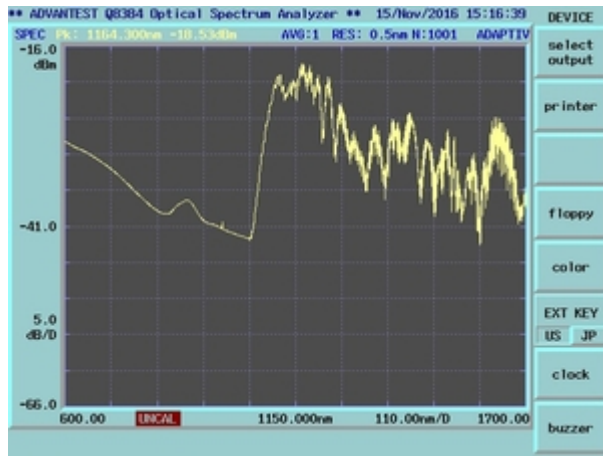


Rf-output-fiber

- measured with optical spectrum analyzer from AG-Morgner
- Resolution: for peak 0.2nm and full range measurement 1.5nm
- pdf-files:

frequency_comb_spectrum_rf-fiber-output.pdf
frequency_comb_spectrum_rf-fiber-output_full_range_.pdf

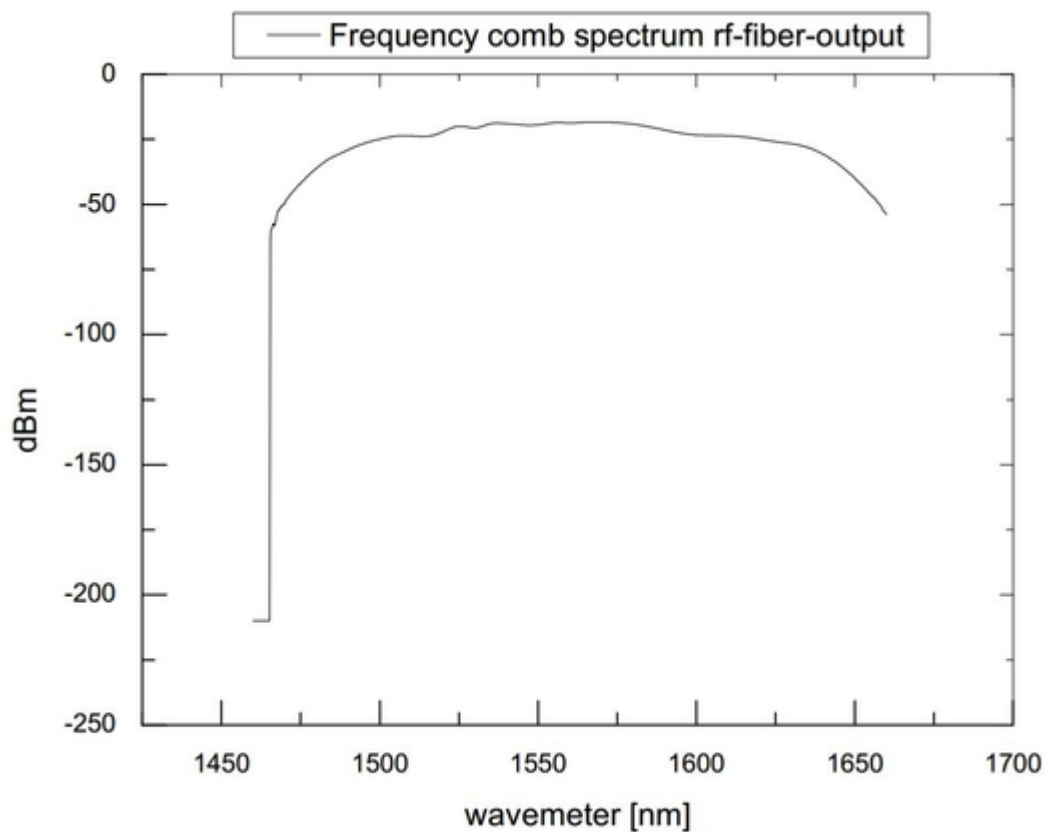


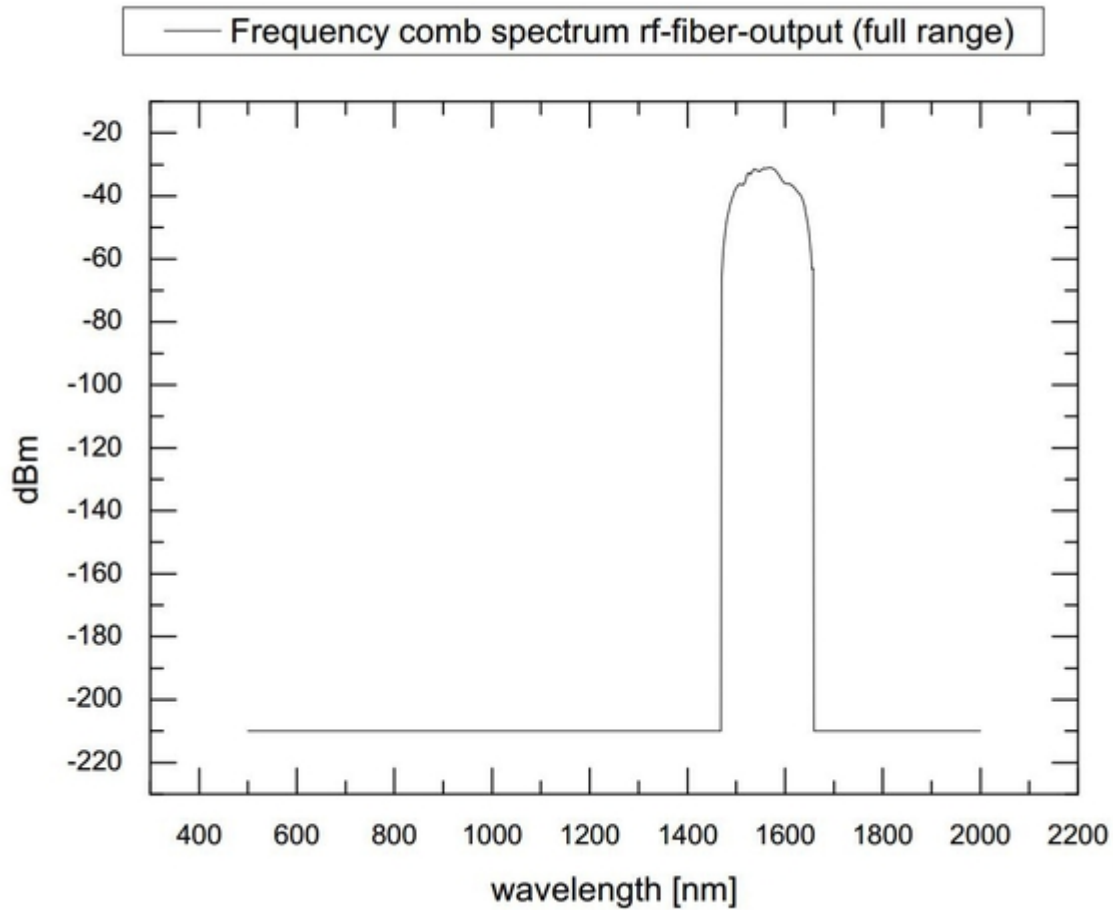



Manuals and Datasheets

- Software for [FiberLaser](#)
- Software for [Fiber Comb Control](#)


Pictures from the inside of the frequency comb

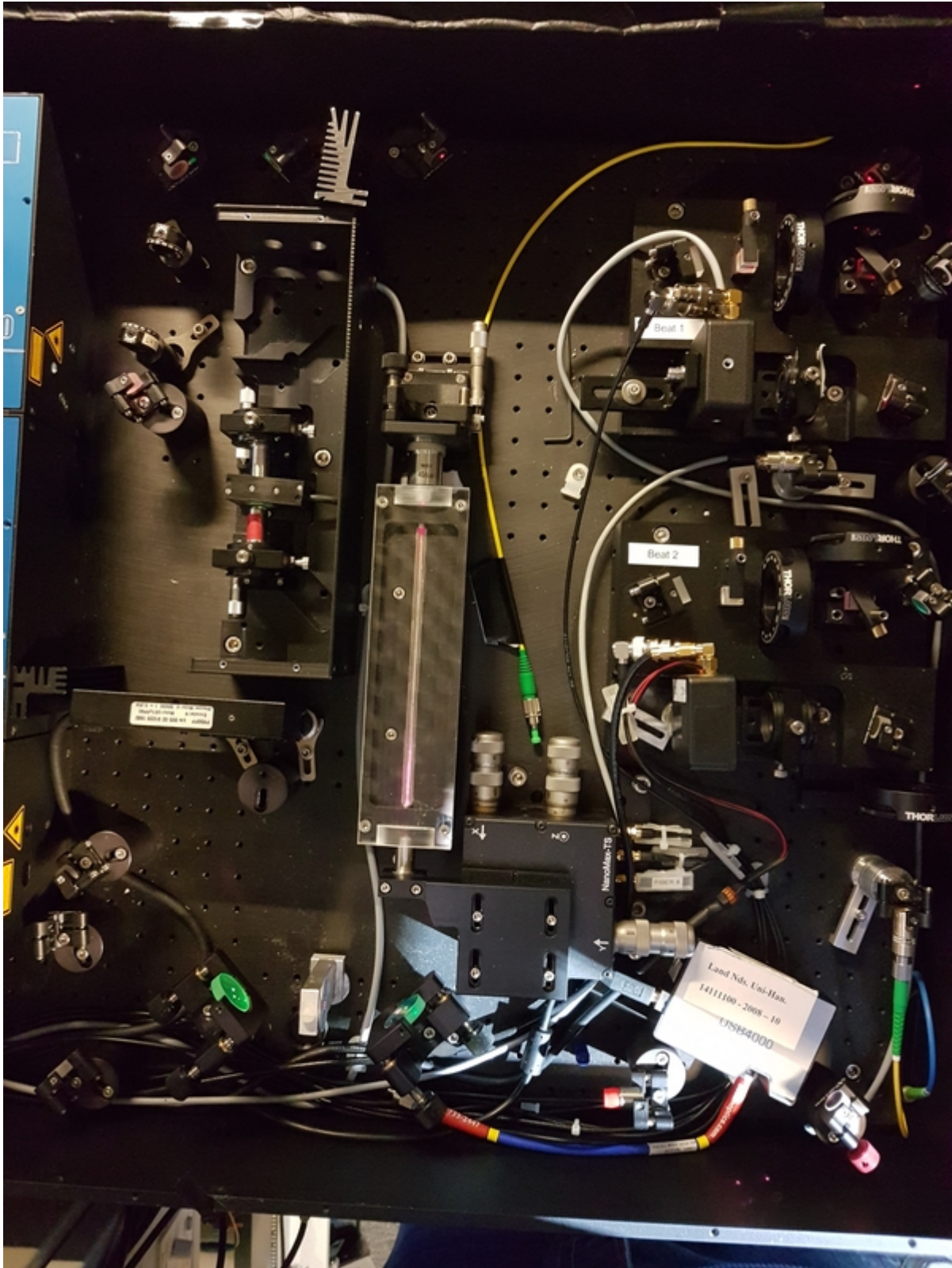


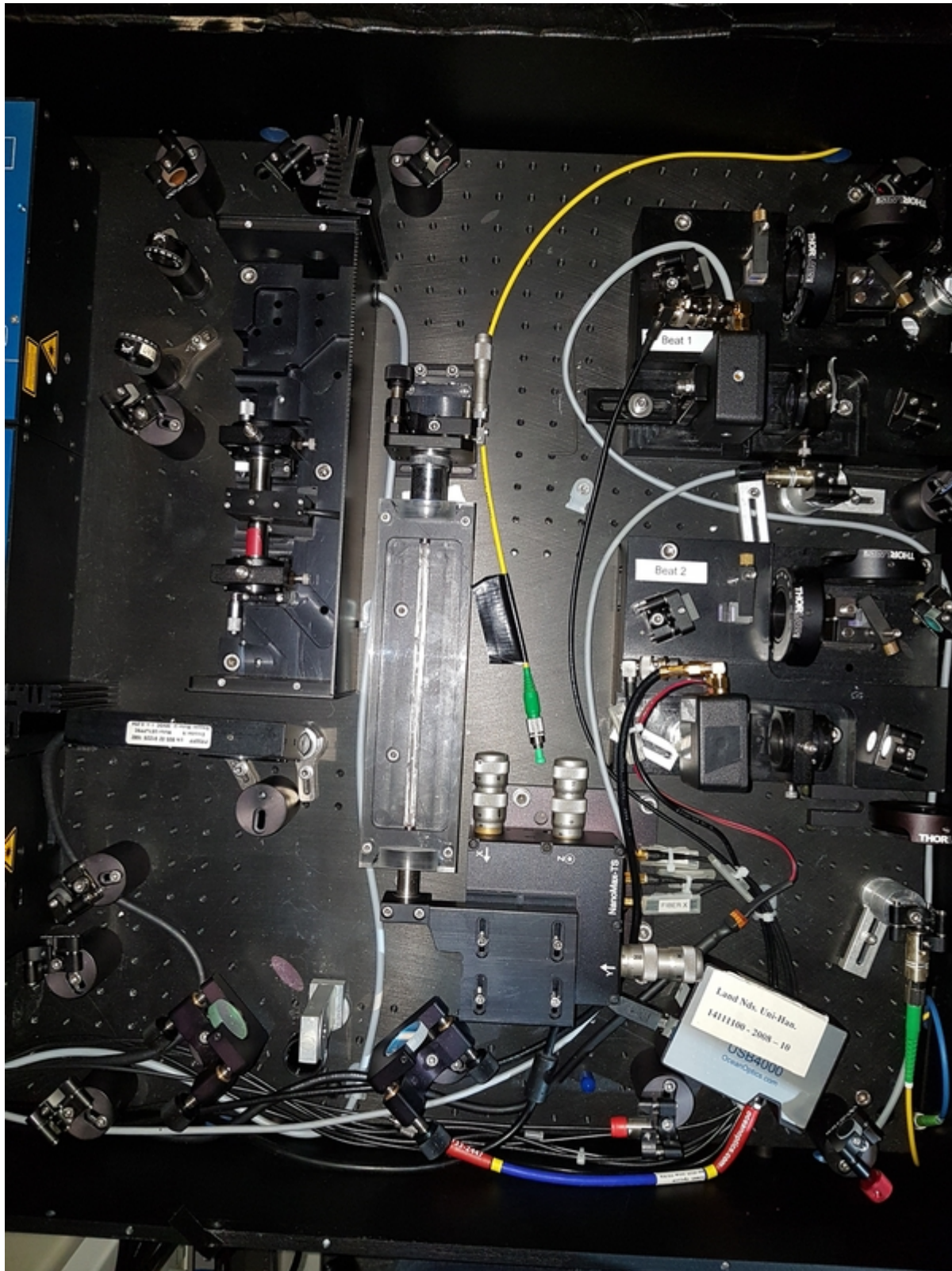


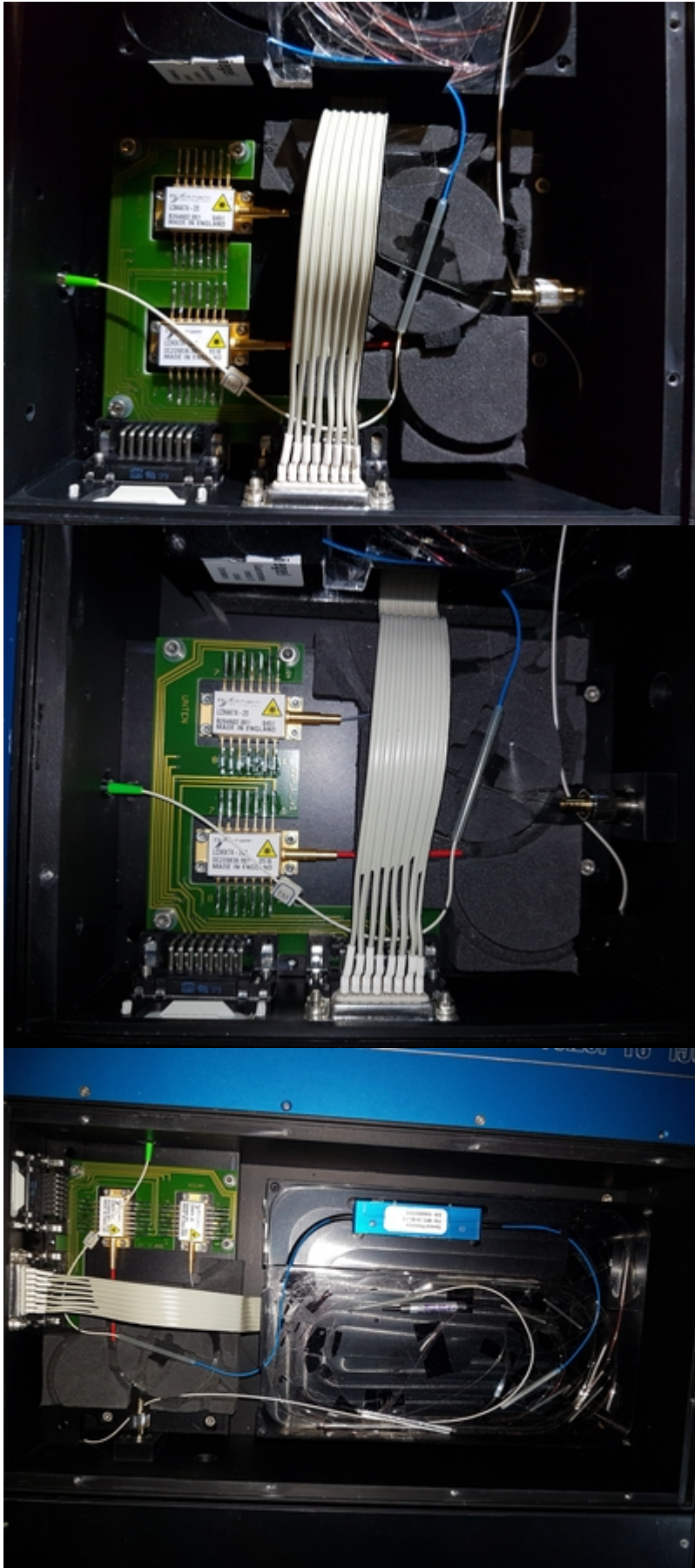





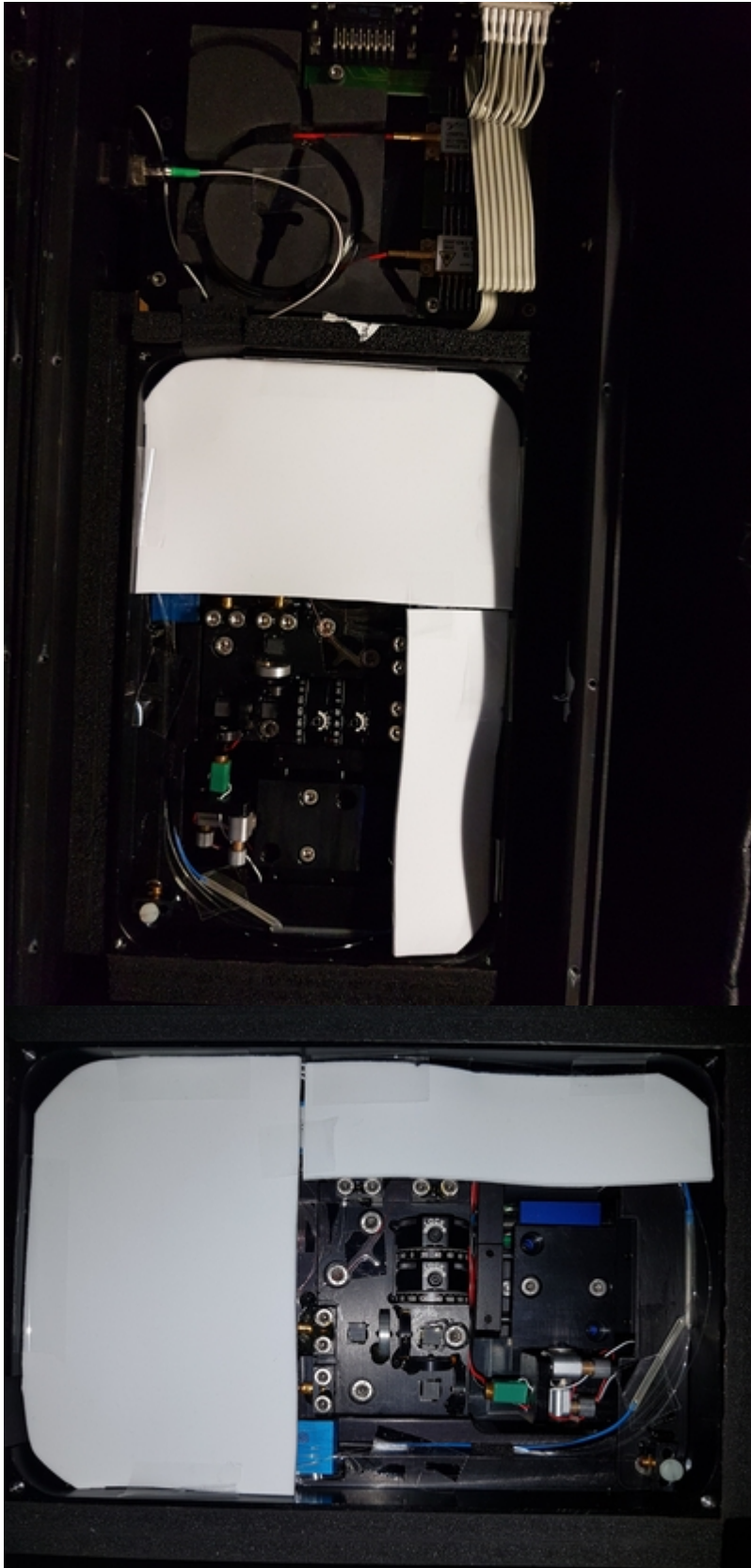



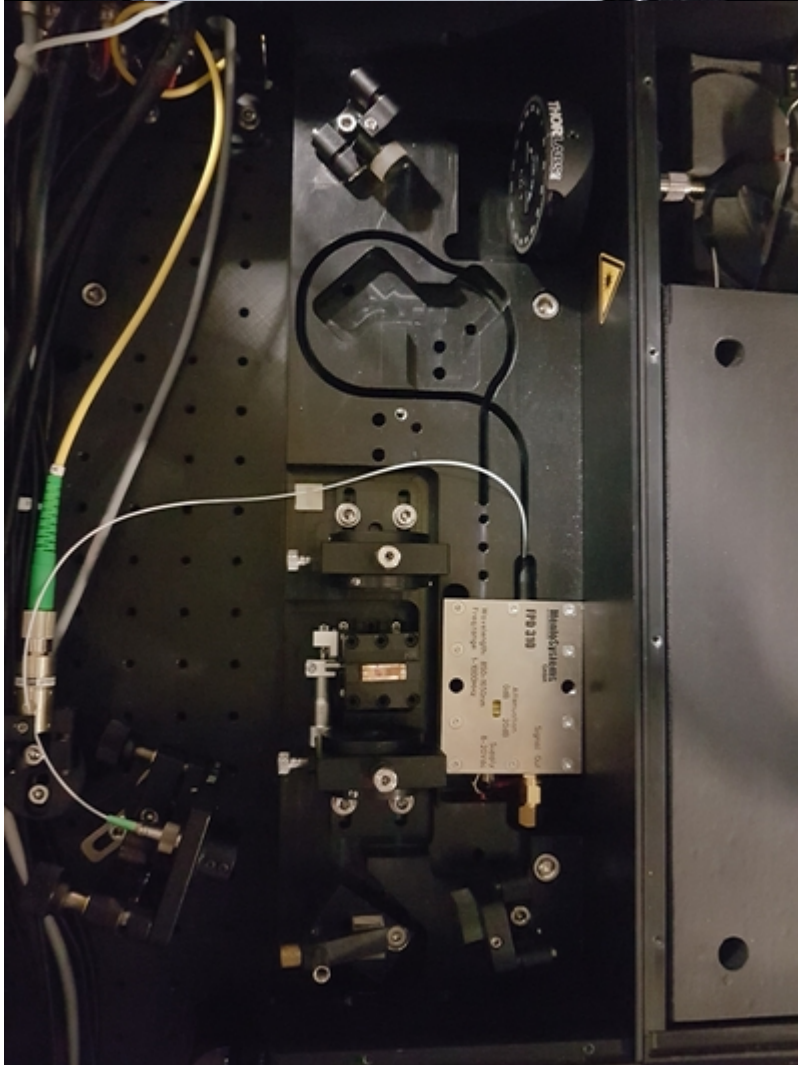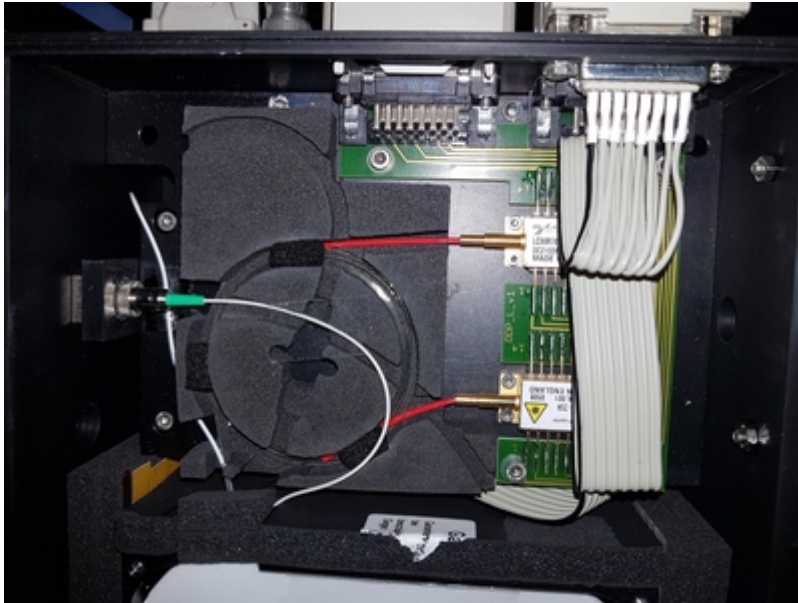

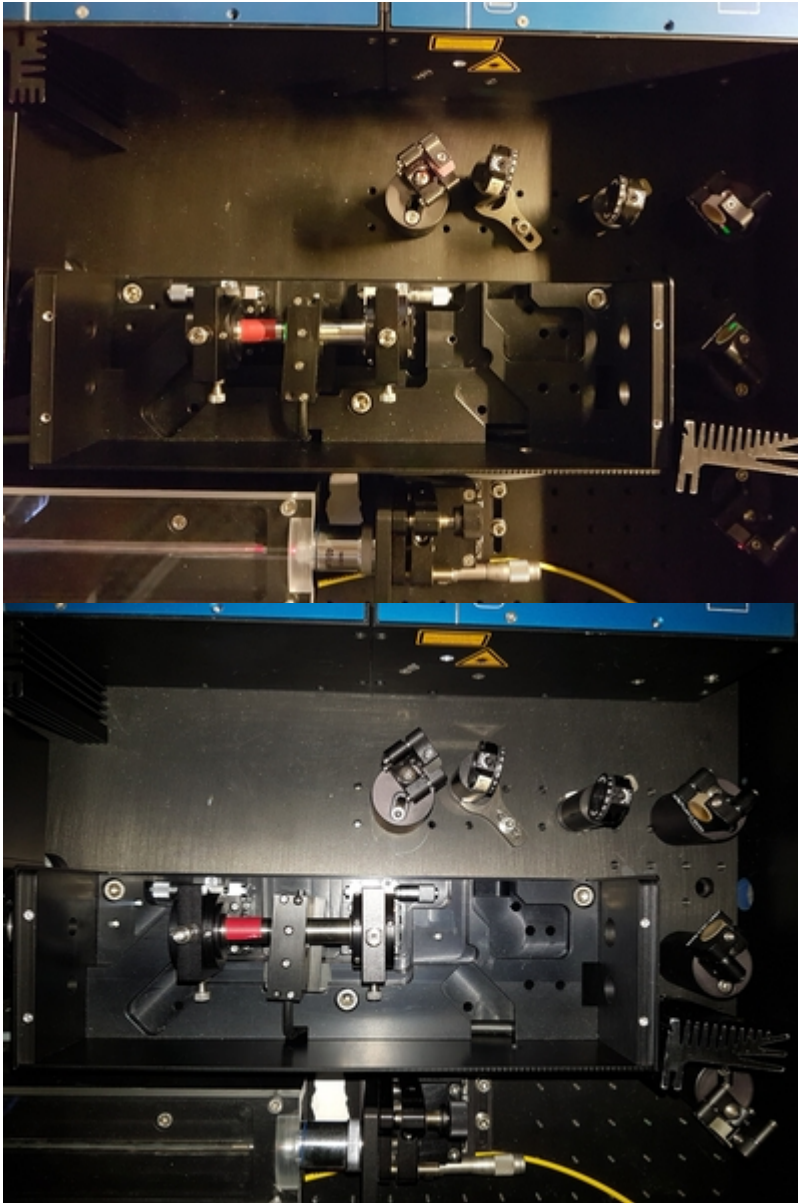

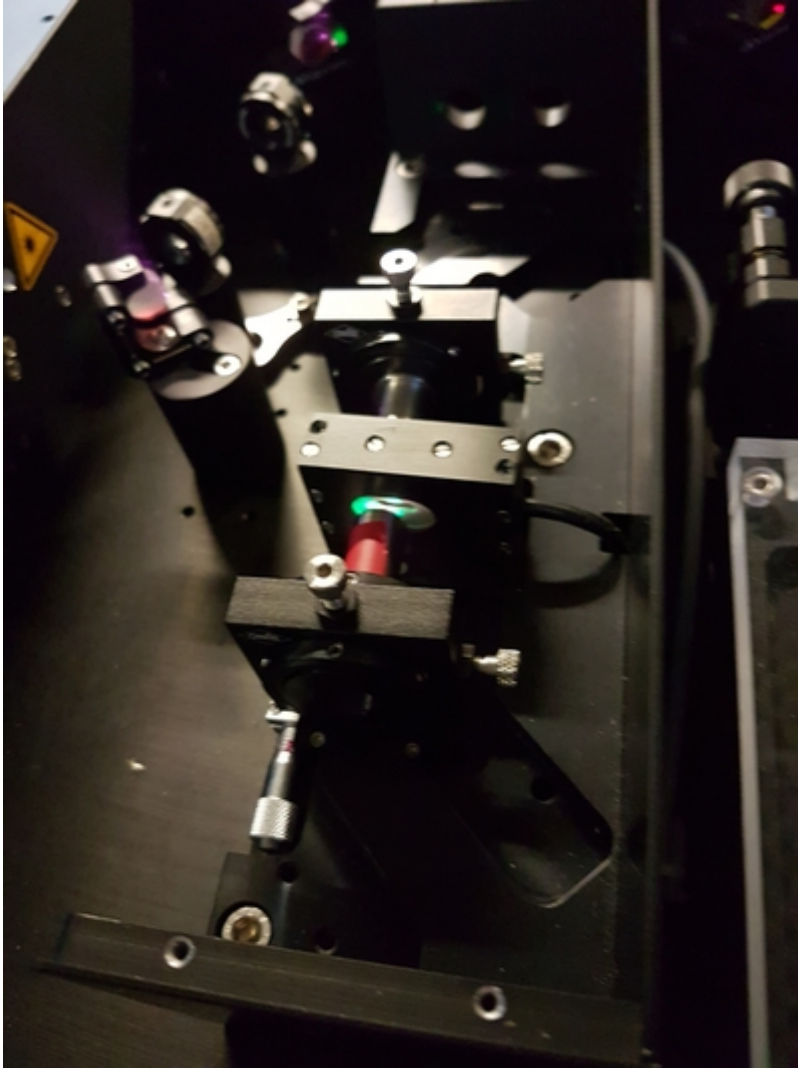

From:
<https://iqwiki.iqo.uni-hannover.de/> - IQwiki

Permanent link:
<https://iqwiki.iqo.uni-hannover.de/doku.php?id=groups:mg:comb&rev=1493295780>

Last update: **2017/04/27 12:23**

